

VACUU•LAN® Mini-Network


from BrandTech® Scientific

The innovative VACUUBRAND® VACUU•LAN® Local Vacuum Network runs up to three applications from a single pump. Not just a simple manifold, the Mini-Network incorporates three chemically-resistant flow-control valves with integrated check valves to minimize interference and cross-contamination. Rather than the typical 100-200mbar (75-150Torr) available from central vacuum systems, a Mini-Network can deliver vacuum as low as 2mbar (1.5Torr) to applications. Powered by an oil-free VACUUBRAND diaphragm pump for quiet reliability, Mini-Networks are a modular way to provide distributed vacuum to a single bench or entire lab with ease.

With the VACUU•LAN® Mini-Network, you can install three VACUU•LAN® ports in a snap! Just attach the assembly to

the wall or your lab frame, connect the pump with vacuum tubing, and you're ready to go! Individual ports are check-valve-controlled to minimize interference among applications. Ports include manual flow control for continuously variable pumping speed adjustment. Coupled with the 100mbar ME1 or ME1C pumps, you have an instant filtration network. Choose the MZ2C NT for higher flow rate filtration applications. For demanding multiple evaporative applications VACUUBRAND offers many pump options including the 7mbar PC101 NT or the 1.5mbar PC201 NT.

Visit www.brandtech.com to learn more about VACUU•LAN® options and use our convenient Vacuum Pump Selection Guide.

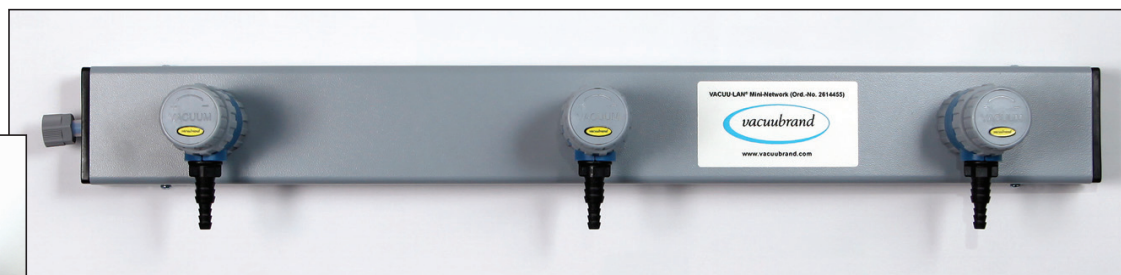

Each VACUU•LAN® Mini-Network includes three manual flow control valves.

ORDERING INFORMATION	Free Air Capacity	Ultimate Vacuum	Cat. No.	2018 List Price
VACUU•LAN® Mini-Network			2614455	\$1,386.52

Suggested Vacuum Supply for Mini-Network

Pump/System

ME1, 100-120V, 50-60Hz, US plug	14lpm	100mbar	721003	\$868.00
ME1C, 100-120V, 50-60Hz, US plug	14lpm	100mbar	721103	1,102.00
MZ2C NT, 100-120V, 50-60Hz, US Plug	38lpm	7mbar	732303	2,978.00
PC101 NT, 100-120V, 60Hz, US plug	38lpm	7mbar	733003	3,791.00
PC201 NT, 100-120V, 60Hz, US plug	63lpm	1.5mbar	737003	7,067.00

ACCESSORIES

Vacuum regulation valve with manometer for ME1			696842	\$182.08
Vacuum regulation valve with manometer for ME1C			696843	365.00
CVC 3000 detect Integrated Vacuum Control System (includes solenoid valve), 100-230V, 50-60Hz, US, CEE, UK, AUS plug			2614120	3,318.00

Product appearance, catalog numbers, prices, specifications, and technical information are subject to change without notice.

BT_MiniNetwork_2018

Lab Rats Trust
BrandTech!



BRANDTECH®
SCIENTIFIC, INC.

Toll Free 888-522-2726
www.brandtech.com