

SUPERIOR PROTECTION FOR LIGHT SENSITIVE MATERIALS


VITLAB® Light Shielding Class A Volumetric Flasks

Protect your light sensitive materials with VITLAB® light shielding volumetric flasks from BrandTech® Scientific. A proprietary pigmented polymethylpentene (PMP) provides a sun protection factor of nearly 20 and is superior to amber glass at virtually all wavelengths, especially for UV light. All with the light weight and superior break resistance of plastic.

- All flasks are individually calibrated to Class A in accordance with DIN EN ISO 1042, imprinted by lot number, and provided with batch certificate.
- Screw caps help protect contents from contamination. Also available with NS stopper
- Can be sterilized by steam (121°C/20 min), formalin, ethanol or β -radiation (25kGy).
- Microwavable.
- Light protection molded in — not a coating, for durability and scratch resistance.

Flasks have excellent resistance to acids, alkalis and alcohols at room temperature, and are suitable for short term exposure to aliphatic hydrocarbons, aldehydes, ketones and esters.

Flasks are not suitable for use with aromatic and halogenated hydrocarbons, ethers, oxidizing acids, oxidizing agents, or other chemicals incompatible with PMP.

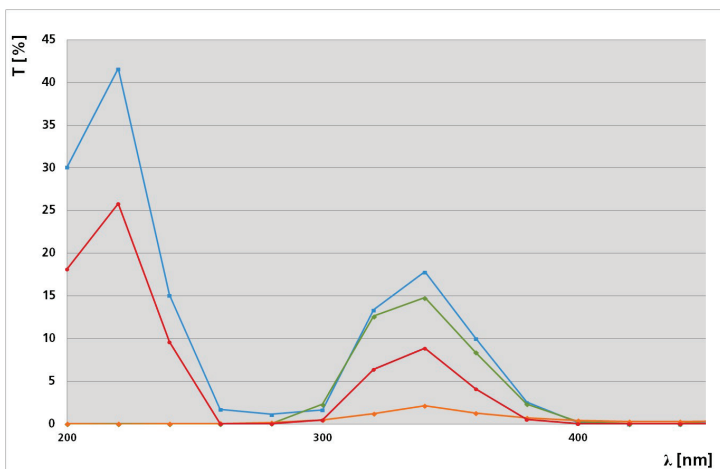

ORDERING INFORMATION

Volume mL	Tolerance +/-mL	Height mm	Thread GL	Qty./Pack	Cat. No.
-----------	-----------------	-----------	-----------	-----------	----------

Class A Light-Shielding Volumetric Flasks, PMP, with Screw Cap, PP, Certified

10	0.04	90	18	2	V670040
25	0.04	115	18	2	V671040
50	0.06	150	18	2	V672040
100	0.10	180	18	2	V673040
250	0.15	235	25	2	V674040
500	0.25	270	25	2	V675040
1000	0.40	310	32	1	V676040

Volume mL	Tolerance +/-mL	Height mm	Stopper NS	Qty./Pack	Cat. No.
-----------	-----------------	-----------	------------	-----------	----------

Class A Light-Shielding Volumetric Flasks, PMP, with NS Stopper, PP, Certified

10	0.04	90	10/19	2	V670950
25	0.04	115	10/19	2	V671950
50	0.06	150	12/21	2	V672950
100	0.10	180	14/23	2	V673950
250	0.15	235	19/26	2	V674950
500	0.25	270	19/26	2	V675950
1000	0.40	310	24/29	1	V676950

Product appearance, catalog numbers, prices, specifications, and technical information are subject to change without notice.

Vitlab_flyer_LightShield-2014



VITLAB

BRANDTECH[®]
SCIENTIFIC, INC.

Toll Free: 888-522-2726
www.brandtech.com