

# Ergonomic Design

## Features of the BRAND® Transferpette® S Pipettes

from BrandTech® Scientific

Transferpette® S single and multichannel pipettes are modern, ergonomic pipettes designed to be easy to operate, while retaining the familiar central piston button design. All controls of the instrument are designed to be operated with small movements of a single hand, reducing the amount of manipulation required compared to other instruments.

### Features a 12.5mm piston stroke length

The Transferpette® S features a 12.5mm piston stroke length – the shortest of all quality pipettes – to reduce thumb movement. Up to 40% shorter than some popular brand

### Volume adjustment brake

A volume adjustment brake can easily be operated by the thumb tip of the pipetting hand, to secure the volume setting from inadvertent drift

### Volume adjustment

Smooth, thumb-tip volume adjustment, which can be secured or released with the brake

### Finger hook


Finger “hook” along with a lightweight and well-balanced design reduces forces required to stabilize the pipette in the hand, resulting in a relaxed hand grip

### Stepped tip ejector

Stepped tip ejector (multi-channel models) dislodges tips two at a time, to spread out the forces required to eject tips

### Upgraded! O-ring design

New O-ring design enhances fit with most major brands of pipette tips. Reduced contact reduces tip mounting forces for a good seal, and likewise reduces ejection forces.

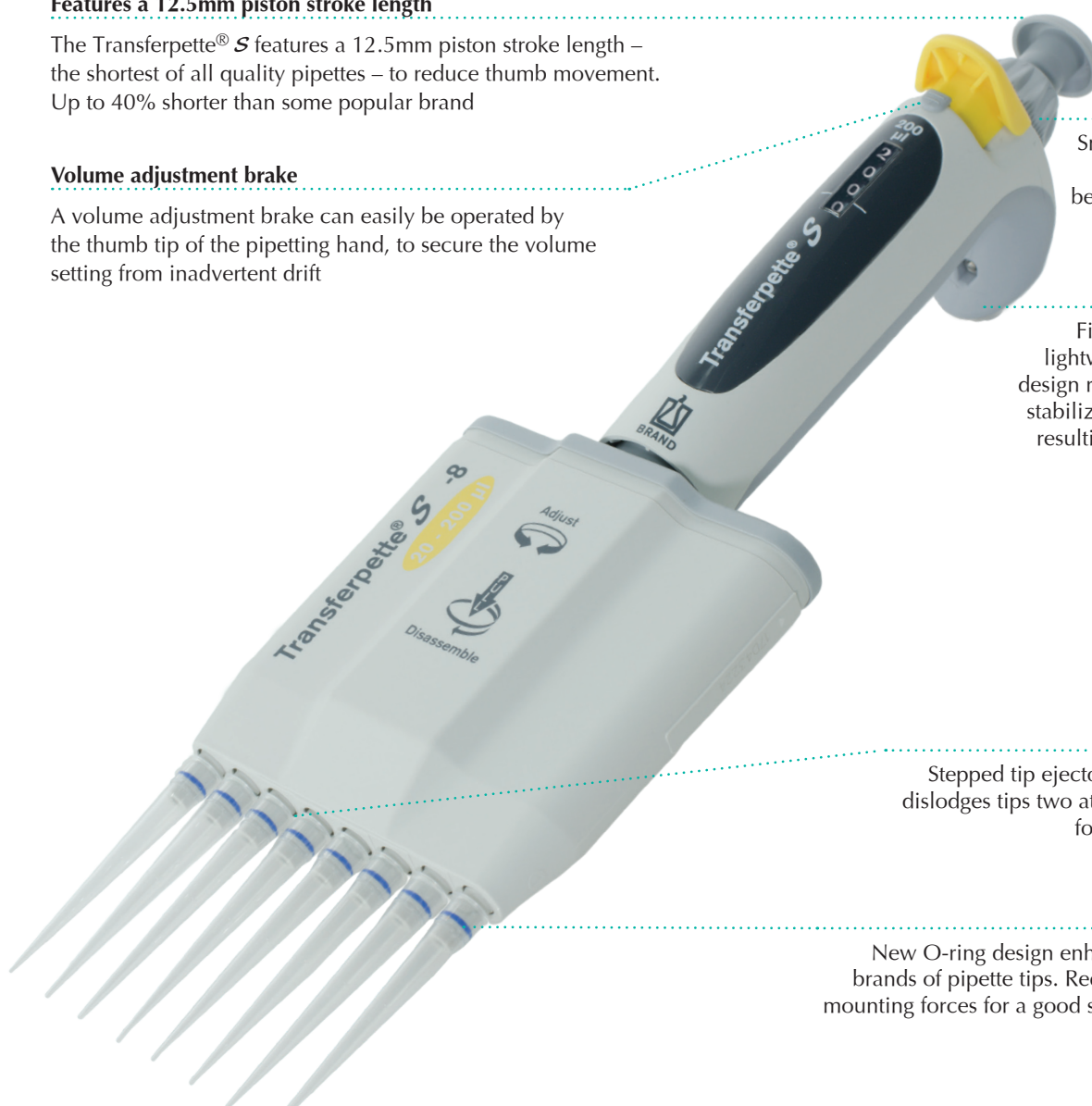

Lab Rats Trust BrandTech®

**BRANDTECH®**  
SCIENTIFIC, INC.

Toll-Free (888) 522-2726

[www.brandtech.com](http://www.brandtech.com)

# Ergonomic Design

## Features of the BRAND<sup>®</sup> Transferpette<sup>®</sup> S Pipettes

from BrandTech<sup>®</sup> Scientific

### Transferpette<sup>®</sup> S Adjustable Volume Single Channel Pipettes

Capacity, $\mu$ L	Model	A* $< \pm$ %	CV* $\leq$ %	Subdivision, $\mu$ L	BrandTech <sup>®</sup> Cat. No.
0.1-1	D-1	2	1.2	0.001	704768
0.1-2.5	D-2.5	1.4	0.7	0.002	704769
0.5-10	D-10	1	0.5	0.01	704770
2-20	D-20	0.8	0.4	0.02	704772
5-50	D-50	0.8	0.3	0.05	704773
10-100	D-100	0.6	0.2	0.1	704774
20-200	D-200	0.6	0.2	0.2	704778
100-1000	D-1000	0.6	0.2	1	704780
500-5000 (0.5-5mL)	D-5000	0.6	0.2	5	704782
1000-10000 (1-10mL)	D-10000	0.6	0.2	10	704784

### Transferpette<sup>®</sup> S Fixed Volume Single Channel Pipettes

Capacity, $\mu$ L	Model	A* $< \pm$ %	CV* $\leq$ %	Subdivision, $\mu$ L	BrandTech <sup>®</sup> Cat. No.
10	F-10	1	0.5	—	704708
20	F-20	0.8	0.4	—	704716
25	F-25	0.8	0.4	—	704720
50	F-50	0.8	0.4	—	704728
100	F-100	0.6	0.2	—	704738
200	F-200	0.6	0.2	—	704744
500	F-500	0.6	0.2	—	704754
1000	F-1000	0.6	0.2	—	704762

Transferpette<sup>®</sup> S Single channel pipettes are supplied with performance certificate, operating manual, shelf mount and three year warranty

### Transferpette<sup>®</sup> S -8/-12 Multichannel Pipettes

Capacity, $\mu$ L	A* $< \pm$ %	CV* $\leq$ %	Subdivision, $\mu$ L	Transferpette <sup>®</sup> S -8 Channel	Transferpette <sup>®</sup> S -12 Channel
0.5-10	1.6	1.0	0.01	703700	703720
5-50	0.8	0.4	0.1	703706	703726
10-100	0.8	0.3	0.2	703708	703728
20-200	0.8	0.3	0.2	703710	703730
30-300	0.6	0.3	0.5	703712	703732

Transferpette<sup>®</sup> S -8/-12 pipettes supplied with performance certificate, operating manual, shelf mount, two tip boxes, reagent reservoir. One-year warranty.



ergo Tfps\_flyer\_2018



Lab Rats Trust BrandTech<sup>®</sup>

**BRANDTECH<sup>®</sup>**  
SCIENTIFIC, INC.

Toll-Free (888) 522-2726

[www.brandtech.com](http://www.brandtech.com)